**M1200**

**LUXURY CLASS**

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## English

Thank you for purchasing this product. Carefully read through these OPERATING INSTRUCTIONS and pay attention to all of the precautions in order to avoid using the product incorrectly.

Keep this manual in a safe place for future reference.

### CONTENTS

<b>SAFETY PRECAUTIONS</b> .....	24
Applicable environment .....	24
User precautions .....	24
Remote control unit precautions.....	24
Conditions of use .....	24
Not to be used under the following conditions .....	25
Assembly and maintenance .....	26
Instructions for the mains cable and plug .....	26
 <b>IDENTIFICATION AND FUNCTION OF EACH PART</b> .....	26
Outer structure .....	26
 <b>PREPARATION BEFORE USE</b> .....	27
Place of installation .....	27
Floor protection .....	27
Earth connection .....	27
Assembly and installation .....	28 – 31
Checking the backrest cushion .....	32
Movement and relocation .....	32
Surroundings .....	32
 <b>MAIN CONTROL FUNCTIONS</b> .....	33
Instructions on the buttons .....	33 - 36
Instructions for the remote control music .....	38
 <b>OPERATION</b> .....	39
Start up .....	39
Before sitting .....	39
Recline and extend system .....	39 - 40
After use .....	40
 <b>FAULT FINDING</b> .....	41
 <b>MAINTENANCE</b> .....	42
 <b>TECHNICAL SPECIFICATIONS</b> .....	43

## ● SAFETY PRECAUTIONS

- Before use, carefully read through the safety instructions and use them correctly.
- Follow the safety instructions to avoid injury or damage to the product.
- For safety purposes to avoid potential damage and misuse of the product, pay strict attention to the instructions strictly marked as follows:

 <b>CAUTION</b>	Improper use can cause serious injury to the user
 <b>WARNING</b>	Improper use can cause serious injury to the user or damage the appliance.

## ● Applicable environments

 <b>WARNING</b>	To avoid electric shock or injury do not use the chair in a bathroom or damp surroundings
	To avoid the risk of fire, do not use the product on or near heating appliances or expose it to direct sun light. Do not use at temperatures above 40°C.
	To avoid the risk of fire, do not use the chair on electrically heated mats or other types of heaters.
	Make sure the chair is placed on a flat surface.
	When changing from cold to hot temperatures, wait one hour to prevent the machine suffering from frost damage.

## ● User precautions

**Caution: Consult your doctor before beginning to use this appliance. This advice is especially important for those suffering from health problems.**

## ● Remote control unit precautions

 <b>WARNING</b>	Do not handle the remote control unit with wet hands
	Do not drop it or immerse it in water
	Do not tug on the connection cable
	Do not place heavy objects on top of the remote control unit
	It is advisable to always keep the remote control unit on its support.

## ❖ Conditions of use

- ▶ To prevent possible injury or incorrect results do not use the machine in conjunction with electric blankets.
- ▶ Do avoid nauseous effects do not use for up to half an hour after eating.
- ▶ do not allow children or household pets to get too near to the machine, backrest, seat or under the leg rest.
- ▶ To prevent damage or injury do not place heavy objects on top of the machine in case they fall off.
- ▶ A massage session lasting longer than 5 minutes in a fixed body position may cause overstimulation and negative effects.
- ▶ The chair is not to be used by more than one person at a time.
- ▶ The chair is designed to be used by people weighing less than 120kg.

	To lower the footrest or backrest, apply moderate pressure and make sure there are no objects obstructing it or children nearby.
	Before use uncover the backrest and make sure that it is in good condition, along with the rest of the component parts. If you detect any damage get in touch with the T.A.S or your local dealer.
	Check that the electric cable is not twisted or caught and do not place heavy objects on it.
	If the mains cable is damaged do not attempt to replace it. This must be replaced by the manufacturer or Technical Assistance Service.
	Do not allow children or disabled people to use the chair without supervision.
	Do not allow children to play with the chair.
	To avoid injury only one person can use the chair at a time.
	Do not use sharp or cutting objects on or near the chair. Otherwise this may damage the chair upholstery.
	Make sure you are not wearing or have sharp or dangerous objects in your pockets, ie cell phones, pocket knives, blades, watches, etc.
	Do not smoke while using the chair.
	Do not use other medical appliances while using the chair.
	Should you experience discomfort or other symptoms during use, stop the chair and consult your doctor.
	After unplugging the mains cable put the electric plug in the place designed to house it.
	When moving the chair make sure that the floors are solid; fragile wooden or carpeted flooring could be damaged.
	Do not use for therapeutical purposes.
	After using the chair, switch it off and unplug it from the mains supply.
	Do not use petrol, caustic cleaning products, bleach or pesticides to clean the chair.
Do not sit on the backrest, armrests or footrest. This could cause injury and damage the chair.	
	Repairs or replacement of parts must always be performed by qualified personnel from the Technical Assistance Service or your local dealer.
	Always unplug the chair from the mains electric supply before storing it away or during maintenance. Make sure your hands are dry.
	In the event of a power cut unplug the chair from the mains supply.

 <b>WARNING</b>	During use remain in an upright position and make sure that the rollers move in unison with your vertebrae to avoid any pain.
	Do not go to sleep while using the machine.
	Do not use if you have consumed alcohol or your body is in any other abnormal state.
	Do not touch the massage rollers directly with your hands, arm joints, knees, chest or abdomen.
	If your legs open out during the massage do not try to force them closed.
	It is advisable to use the chair 20 minutes at a time, and for maximum of 5 minutes in each position.
	To avoid injury only one person can use the chair at a time.

### ❖ Not to be used under the following conditions

If any of the following situations occur, switch the chair off and unplug the mains cable:

- ▶ In case of malfunction, contact your local dealer or technical service.
- ▶ Accidentally splashing water on the chair, the control unit or auxiliary remote control unit.
- ▶ Any of the internal components or accessories become exposed.
- ▶ If the user experiences any problems during use.
- ▶ Any faulty functions during use.
- ▶ The chair suddenly stops.
- ▶ Thunderstorms.

## ❖ Assembly and maintenance

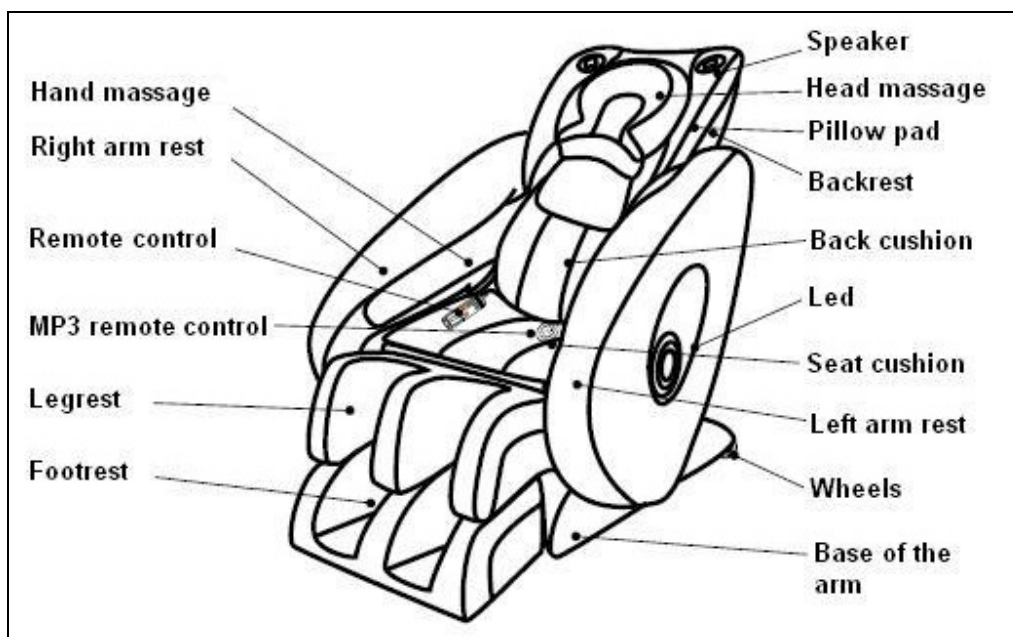
- ▶ Do not remove any of the leather covers to gain access to internal components. Danger of electric shock.
- ▶ In case of problems, contact an authorised technical service. Do not attempt to dismantle the machine yourself.

## ❖ Instructions for the mains cable and plug

- ▶ Check the rated voltage
- ▶ Do not unplug or plug the machine in with wet hands as this may cause electric discharge or damage to the product.
- ▶ Do not tug on the cable to unplug the chair, remove by the plug.
- ▶ The selected timer functions may cause the chair to switch off automatically. Make sure to switch the chair off after use.
- ▶ Do not place heavy objects on the mains cable.
- ▶ If the cable suffers damage, stop using the product, switch the chair off and unplug it from the wall socket. Get in touch with your technical assistance service to have it repaired or replaced.

## ● IDENTIFICATION AND FUNCTION OF EACH PART

### ● Outer structure

## ● PREPARATION BEFORE USE

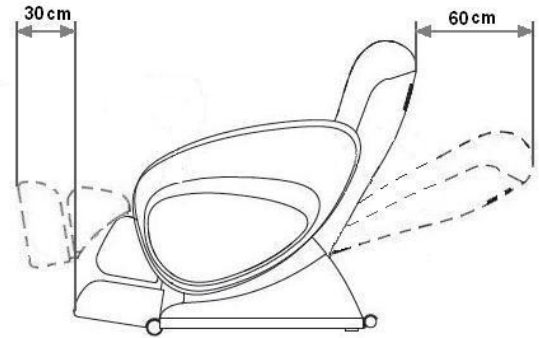
### ● Place of installation

- Safety space around the chair.

At least 60 cm for the backrest between its highest and lowest point.

At least 30cm for the footrest between its highest and lowest point.

To avoid electrical interference the chair must be at least 1m away from television sets and radios.


### ❖ Floor protection

This chair may damage the floor, so it is advisable to place a mat, for example, under the chair to protect the floor and reduce the noise.

Mat

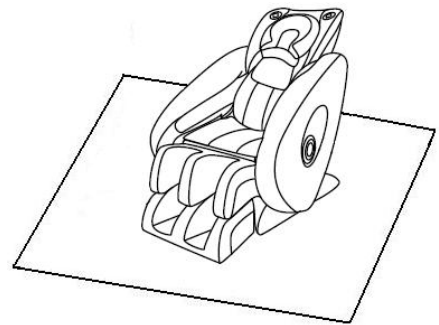Floor

### ❖ Earth connection instructions:

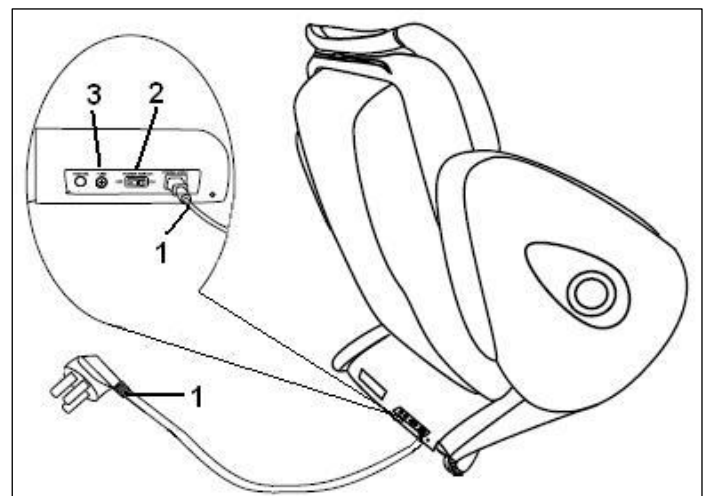
This is a Class I electrical appliance.

Be sure to check the earth connection to prevent electrical leakage or any accidents during use.

An inadequate earth connection could result in electric shock.

Do not modify the electric plug supplied with the massage chair. If the plug is not suitable for the wall socket, contact an electrician to install a properly earthed plug on massage chair.


Insert the plug (1) fully into a suitable socket.

## ❖ Assembly and installation

- Unpack the chair and footrest separately, this will require a minimum of four people. Check the parts list to make sure nothing is missing.

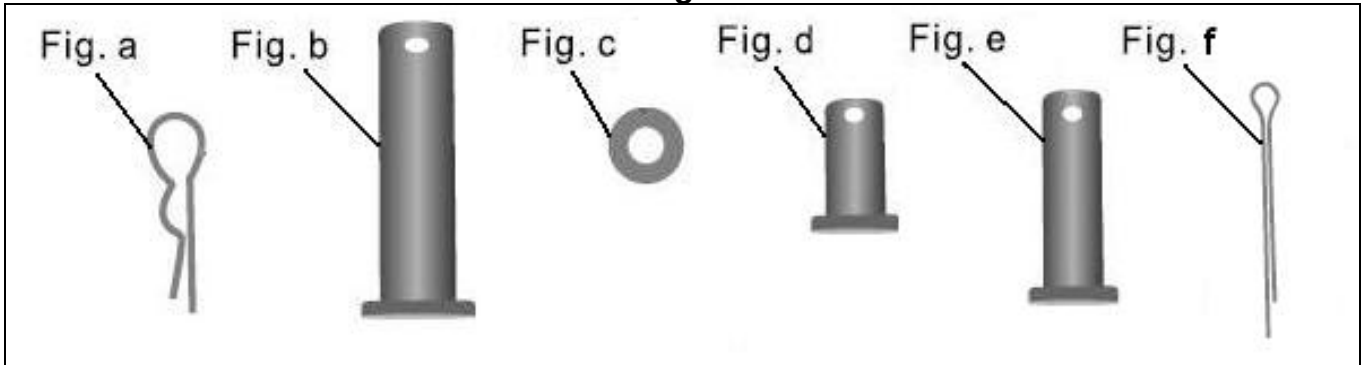**Fig.1**

**Bag 1.-**

**Fig.a.-** Cotter pin.

**Fig.b.-** Arm motor fastening shaft.

**Bag 2.-**

**Fig.c.-** Flat washer

**Fig.d.-** Legrest fastening shaft.

**Fig.e.-** Legrest joint fastening shaft.

**Fig.f.-** Split pin.

► **Step 1:**


After having unpacked the chair, insert the plug (1) into the mains socket (220v), Fig.2.

Set the ON/OFF switch to position (1).

Press the start  button on the remote control unit.

Press the raise backrest button and bring the backrest up to its maximum position.

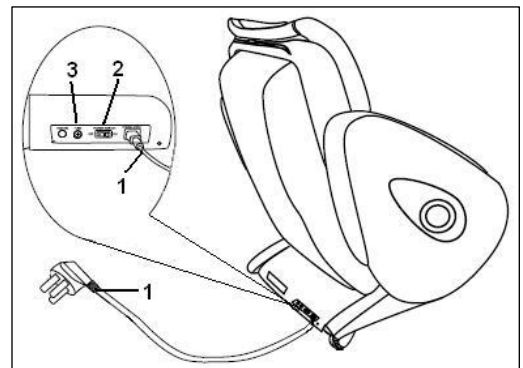
**Fig.2**

► **Step 2:**

Take the screw (M) off the chair's seat cover, as shown in Fig.3.

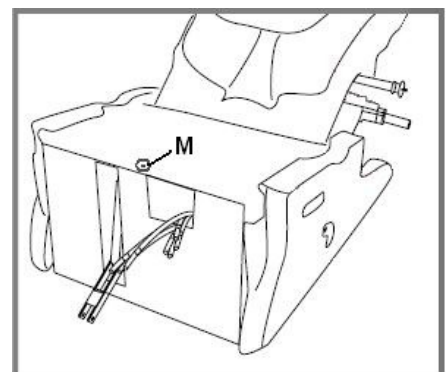
**Fig.3**


► **Step 3:**  
**Fitting the left armrest.-**

Remove the screw and washer (Fig.K) from the top shaft protruding from the chair, Fig.4.

Press the **Arm open** button on the remote control, Fig.5, to bring the arm adjustment shaft out to its maximum position, Fig.L, see Fig.4.

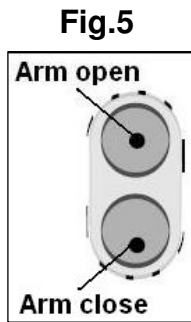
**Fig.4**

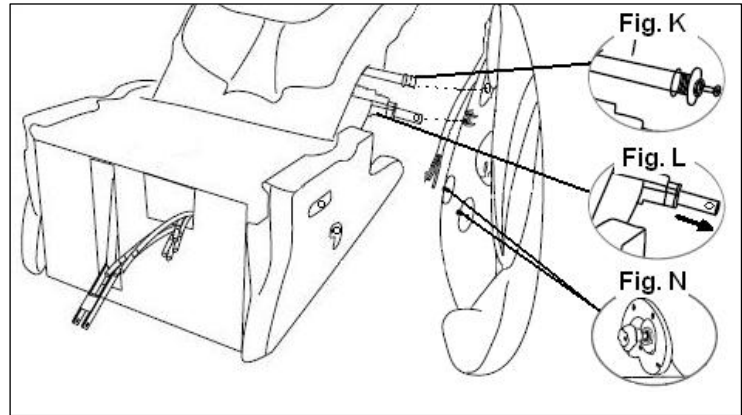
Take the knobs and outer flat washers Fig.N off the armrest, Fig.4.

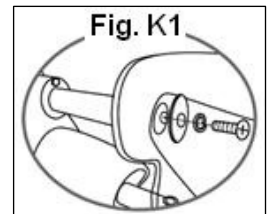Fit the lefthand armrest by slotting it onto the shaft on the chair, Fig.K1, then replace the screw and the washers removed previously and tighten.

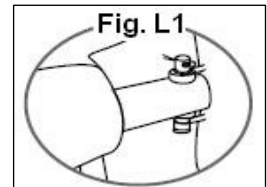
Next, bring the armrest back in toward the chair, insert the coupling pin in through the motor shaft, Fig.L1.

Fit the arm motor fastening shaft Fig.b and insert the cotter pin Fig.a.

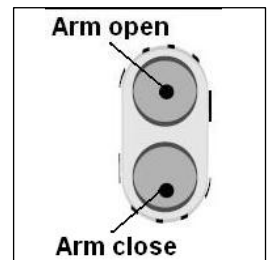Press the **Arm close** button on the remote control to bring the arm adjustment shaft in and be able to insert the armrest shafts, Fig.N, into the holes on the chair, Fig.4. Replace the flat washers and knobs removed previously and tighten.

► **Step 4:**  
**Fitting the right armrest.-**

To fit the righthand armrest, follow the same assembly procedure as done with the lefthand armrest.

**IMPORTANT:**


Press the Start/Stop button on the remote control, press the ON/OFF switch, Fig.2 and remove the plug from the wall socket to proceed with the chair assembly.

► **Step 5:**  
**Fitting the air pipes and the electrical terminals for the armrests.**

**Fig.5**

After having fitted the left and right armrests, connect the lefthand armrest air terminals, marked with the letters (F- G- H- I) and the electrical terminal (AP), Fig. 5.

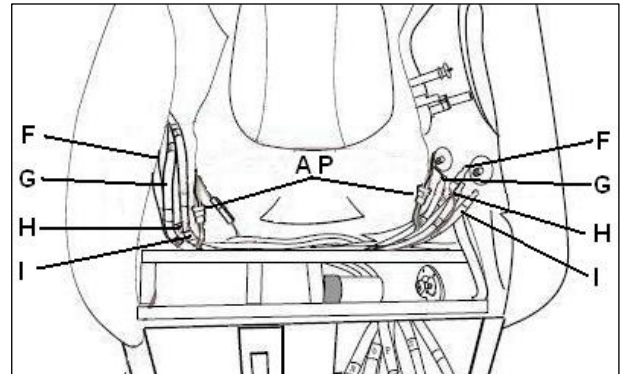Next, connect the righthand armrest air terminals, marked with the letters (F- G- H- I) and the electrical terminal (AP), Fig.5.


► **Step 6:**  
**Fitting the leg support.-**

First off, take the two sets of three and four air pipes, marked (O- P- N) and (Q- R- S- T), out from inside the chair, Fig.6.

Fit the seat cover which was removed in step 2, Fig.3, using the screw that was removed previously.

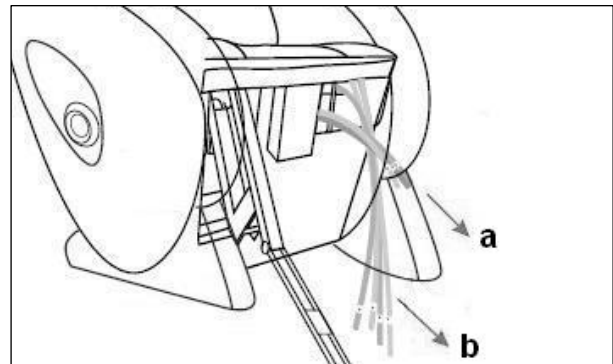
**Fig.6**

Next, lift the swivel joint in the direction of the arrow, Fig.7 and fit the legrest fastening shaft Fig.d.


Insert the split pin Fig.f and bend the ends over as shown in Fig.7.

**Fig.7**

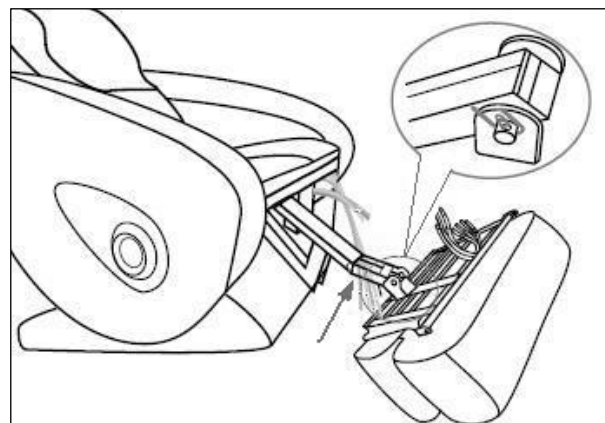
**Fig.8**

Connect the two sets of air terminals, as shown in Fig.8.

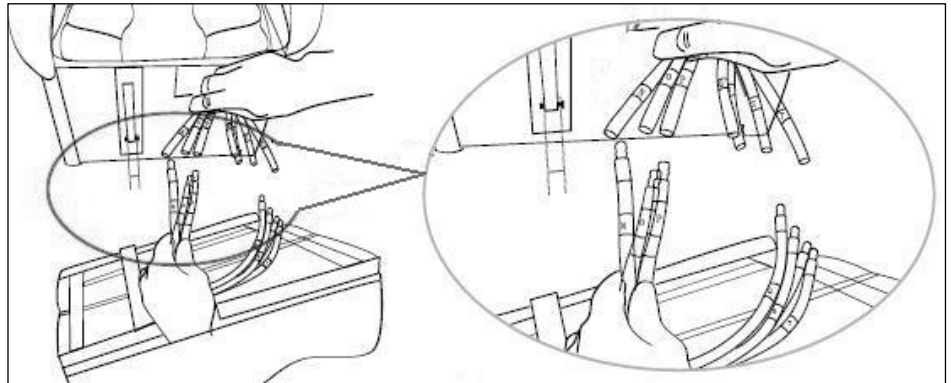**Fig.9**

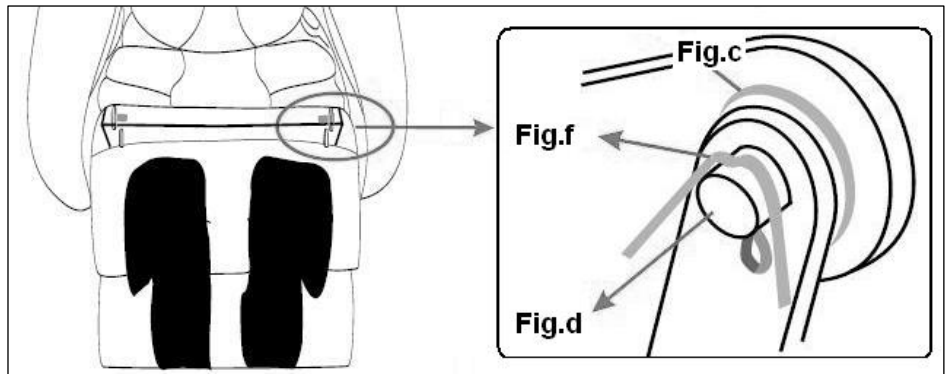Next, rest the footrests on the two wings of the seat, Fig.9. Insert the legrest fastening shafts Fig.d along with the washer Fig.c and insert the split pin Fig.f, bend the ends of the split pin over as shown in Fig.9.

► **Step 7: Attaching the cushion pad:**

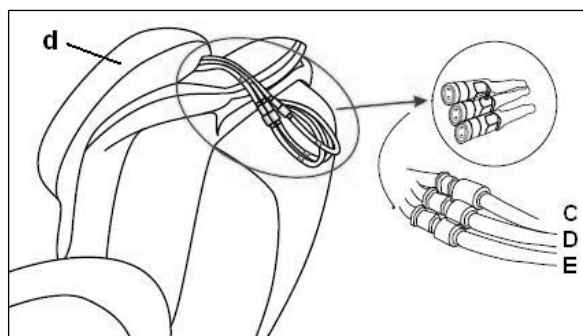
**Fig.10**

Position the cushion pad (c) and, after connecting the right and left audio terminals, close the zip as shown in Fig.10.

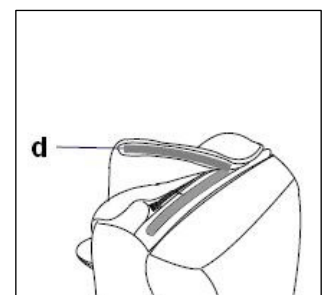
**Fig.11**

Next, position the head massager (d) at the top of the backrest, connect the quick-fit air connectors (C-D-E) Fig.11 and use the Velcro to hold it in position, as shown in Fig.12.


**Fig.12**

## ❖ Checking the cushion pad:

- Lift the cushioning mat Fig.13 to inspect the backrest cushion.
- If you detect any holes or cracks in the cushion, stop using the chair and get in touch with your technical assistance service for immediate repair.

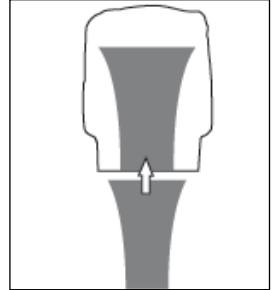**Fig.13**

## ❖ Movement and relocation.-

The help of a second person is recommended when moving the chair.

1. Lift the legrest.

2. Next, use both hands to tilt the front of the chair back on the two rear wheels, approximately 30 degrees, as shown in the Fig.14.


### **Caution:**

1) Move the unit over smooth surfaces; the chair is heavy so be very careful when moving it.

Otherwise you could suffer an injury.

2) Do not move the chair with anyone or anything on it, as it may topple over.

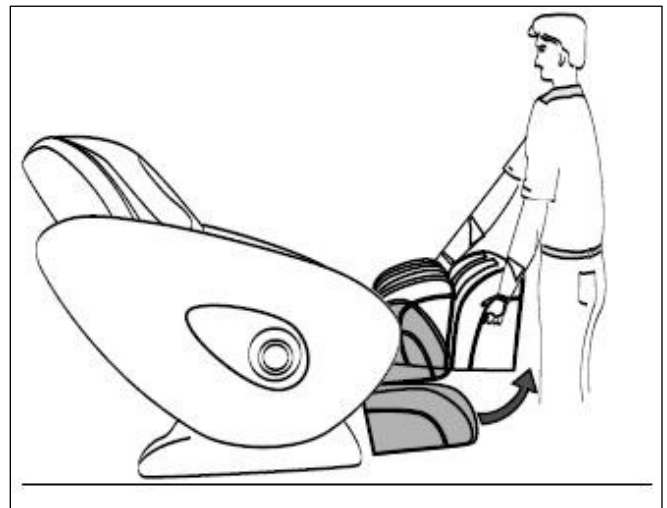
**Fig.14**

## ❖ Surroundings:

- Do not place the chair in the bathroom or anywhere damp, as this could cause electric shocks Fig. 15.
- Do not place the chair near heaters or heating appliances as this could cause a fire or damage the leather upholstery.

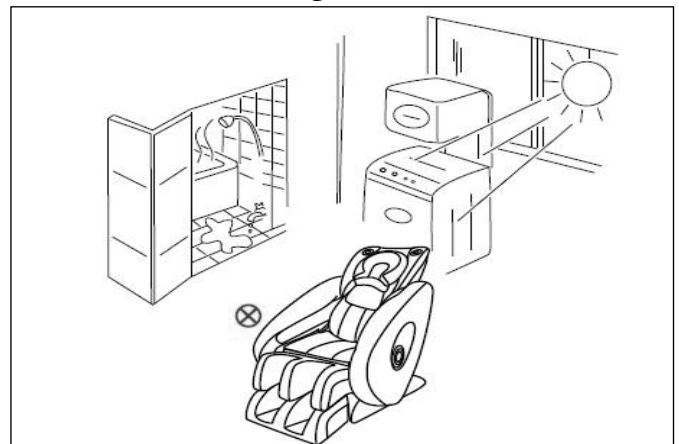**Fig.15**

## ❖ CONTROL PANEL FUNCTIONS.

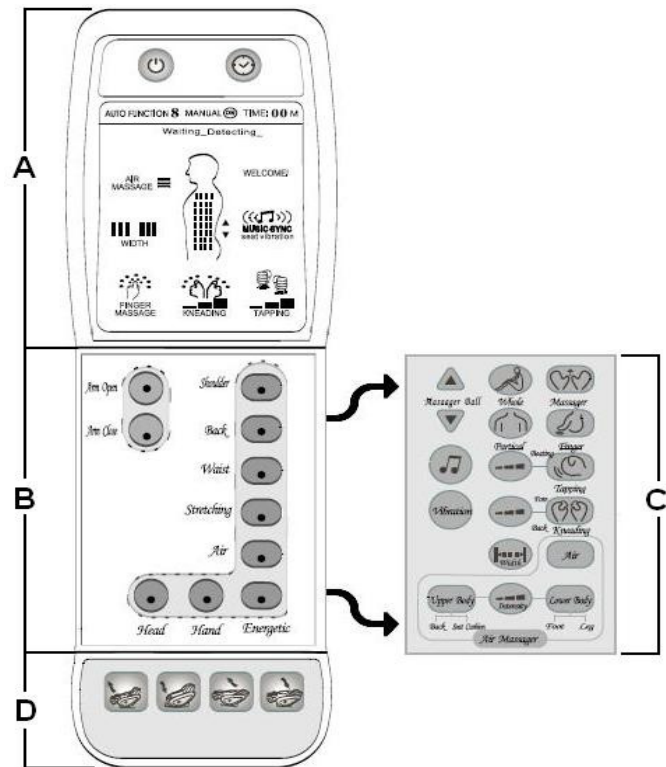
The control panel comprises of four sections with keypads.

**A)** Top keypad and function display.

**B)** Outer mid keypad, automatic and manual programs.

**C)** Inner mid keypad, massage settings and modes.

**D)** Bottom keypad for chair movements.

### ● Key/button functions.-

**ON/OFF.**

Press once to switch on the unit.  
Press it again to switch it off.

**TIME.-** Time settings

The default setting is 15 minutes. Each time this key is pressed it adds 5 more minutes, up to a maximum of 30 minutes.

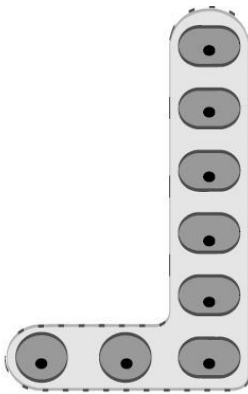
**ARM OPEN.-** open armrest

**ARM CLOSE.-** close armrest

Press and hold this key to open the armrest out. The armrest will stop when the button is released.

Press and hold this key to close the armrest in. The armrest will stop when the button is released.

**SHOULDER.-** Shoulder program.

Automatic program for neck and shoulders. After selecting this option the unit will mainly massage the neck and shoulders.

**BACK.-** Back program

Fully automated back massage program. After selecting this option the unit will mainly massage the back to relax it.

**WAIST.-** Waist program

Fully automated waist massage program, after selecting this option the unit will mainly massage the waist and hip area to relax them.

**STRETCHING.-** Stretching program.

Fully automated stretching program to relax bones and muscles.

**AIR.-** Air pressure program.

Fully automated air pressure program. After selecting this option the unit will use air pressure to massage the feet, hips and head, as well as changes in pressure to stimulate the muscles and eliminate fatigue.

**ENERGETIC.-** Energy program

Fully automated vitality program. After selecting this option the unit massages to eliminate fatigue in the body.

**HAND.-** Manual hand program

Manual hand program. After selecting this option the unit will use air pressure to massage and relax the hands.

**HAND.-** Manual head program

Manual head program. After selecting this option the unit will use air pressure to massage and relax the head.

**WHOLE.-** Full body massage

Automatic full body massage, systematic massage of the entire body.

**PARTIAL.-** Random massage setting

To activate this function press the button once.

**MASSAGER BALL.-** Raise and lower the rollers.

After selecting the random massage function, you can then select the desired position by using the “▲” key or the “▼” key based on the position of the massager ball. Once the position is established the massage motion will be circular.

Use the “▲” or “▼” keys again to change the position.

**MASSAGER.-** Automatic massage.

Press the key once to start the auto-massage function. Press it again to stop it.

**FINGER.-** Finger massage.

Press the key once to activate this function. Press it again to switch it off.

**MUSIC SYNCHRONIZATION.-**

Press the key once to access the music synchronisation function and be massaged in time with the rhythm selected. Press it again to switch this function off.

**TAPPING.-**

Press the key once to activate this function. Press it again to switch this function off.

**TAPPING STRENGTH.-** Tapping force settings

When this function is activated, medium speed is selected. Press the key once and it will remain at medium speed, twice faster and a third press slower. The speed changes in medium-fast-slow cycle.

**VIBRATION.-** Seat vibration

Press the key once to switch on seat vibration. Press it again to switch this function off.

**KNEADING.-** Forehand and backhand massage.

Press the key once to activate forehand massage, twice for backhand massage and three times to switch this function off.

**FORE / BACK.-** Massage setting

Press this key once for medium speed, twice for fast and three times for slow speed.

**WIDTH.-** Roller width setting

**Note:** This key can only be used to adjust the settings when the finger or tapping massage modes are active (wide, medium, narrow).

## Air pressure key set.-

### **AIR.-** Air pressure

Press this key to massage with the following keys:

If you prefer not to use the mechanical massage, press the up “▲” or down “▼” key to disable the mechanical massage function.

### **UPPER BODY.-** Air pressure on the top part of the body

Press this key and massage will be applied to the top part of the body, press it again to massage the back, press it a third time to shift the massage to the cushion area and a fourth time to switch the function off.

### **LOWER BODY.-** Air pressure on the lower part of the body

Press this key once to begin a foot and leg massage, twice for the feet only and three times for legs only. Press the key a fourth time to switch this function off.

### **INTENSITY.-** Air pressure intensity settings

This displays the intensity of the automatic air pressure program and the air pressure setting.

### **LIFTING FOOT SUPPORT.-**

This button electrically raises the footrest. Press and hold the button to activate the function. The footrest will stop when the button is released.

### **LOWERING FOOT SUPPORT.-**

This button electrically lowers the footrest. Press and hold the button to activate the function. The footrest will stop when the button is released.

### **LIFTING BACKERS.-** Raise backrest

This button electrically raises the backrest. Press and hold the button to start the function, which will also lower the footrest. The footrest will stop when the button is released.

### **LOWERING BACKERS.-** Lower backrest

This button electrically lowers the backrest. Press and hold the button to start the function, which will also raise the footrest. The footrest will stop when the button is released.

## ❖ Principles of massage.-

Gravity and poor posture when standing or sitting cause back problems in 90% of the population. Failure to follow proper treatment in time, may lead to other problems.

	<p><b>Massage:</b> This is the massage chair's essential function. As with traction treatment in hospital, the massage chair will stretch the spine to reduce pressure on the nerves and restore normal operation to body organs.</p>
	<p><b>Finger pressure:</b> Because human muscles are tense for long periods, this reduces blood circulation and causes muscle pain. The finger pressure function is like a soothing massage performed by a therapist. The adjustable force only works on the muscle under the skin to achieve muscle relaxation, improve blood circulation, eliminate pain and relieve fatigue.</p>
	<p><b>Kneading:</b> Long working hours without physical activity causes pain in the back and neck muscles. The massage function eliminates back and neck muscle pain and improves blood circulation.</p>
	<p><b>Tapping and kneading:</b> When our body is exhausted and we suffer severe muscle pain, tapping and kneading are applied at the same time. It is a kneading function with a sensation of tapping. This function provides deep stimulation to acupuncture points for painful muscles and eliminates fatigue.</p>
	<p><b>Air pressure:</b> Circulation problems can cause heart disease, heart attacks and angiosclerosis. We should walk a minimum of 10,000 steps a day to achieve normal blood circulation. Frequent use of air pressure massage with air discharge allows normal blood flow to the heart, and reduces the chances of heart disease.</p>


## ❖ Music control panel

- Built into the unit are speakers, located at the top of the backrest, Fig.16, so you can enjoy your favourite music.

Plug the USB cable into the USB socket.

After 5-10 seconds you can use the remote control unit to listen to the music.

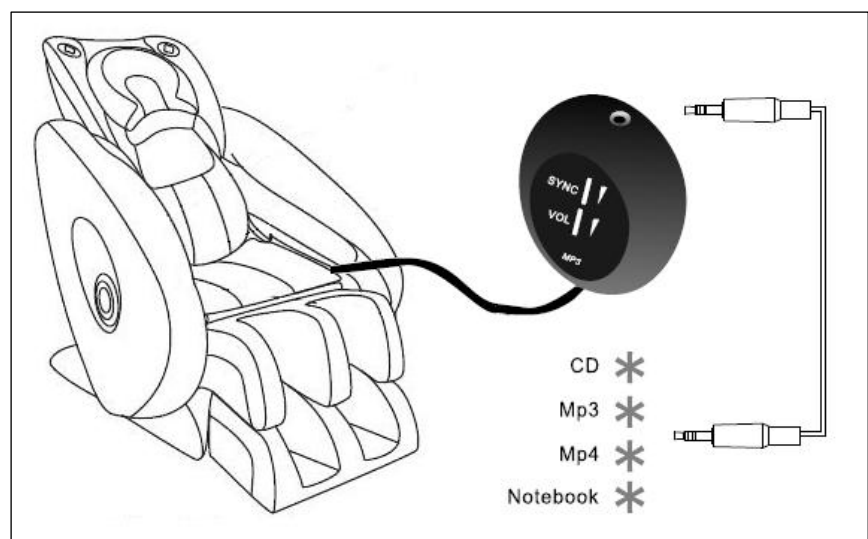Fig.16


## ❖ Instructions for the music control panel

### • Keys

1).- **SYNC.**- Music channel selection

2).- **VOL.**- volume control

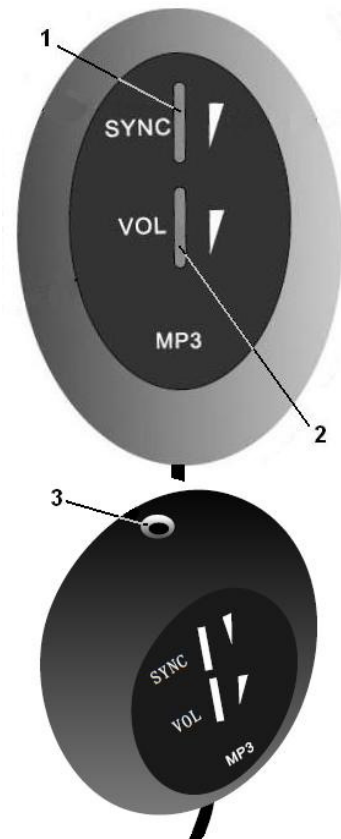3).- Audio input socket

**NOTE:** The volume setting may change the effect of the rhythm massage mode.

## ❖ Music synchronisation.-

The M-1200 massage chair has integrated and adopted ergonomics and AV synchronisation to achieve maximum relaxation for body and mind. You can enjoy both sensory and listening pleasure. The chair allows you to select different massage modes while listening to your favourite music

Feeling of comfort	Sensation of massage + music + synchronisation + video and audio
	Sensation of massage + music + musical synchronisation
	Massage + music
	Just massage

The customer can enjoy musical synchronization by connecting the chair to their laptop, CD, VCD (audio line out) and Mp3 player.

1. Place the MP3 player into the pocket on the righthand armrest. Insert the audio cable into the audio socket.  
Laptop computer, CD or VCD can also be used.
2. Use "Music sync" to select favourite tunes.
3. To adjust the volume control, move the lightly toward "High" or "Low".

**Tip:** When using the "Music Sync" button, adjust the position slowly to achieve optimal synchronization.

## ❖ OPERATION


### ❖ Start up

- Plug the mains cable (1) into the mains supply, Fig.17.
- Set the switch (2) to the ON position.

### ❖ Before sitting

- Before sitting, make sure that all of the components are attached correctly to the chair frame. Make sure that all of the components are in the default position. If not, switch the chair on again to reset all of the positions.

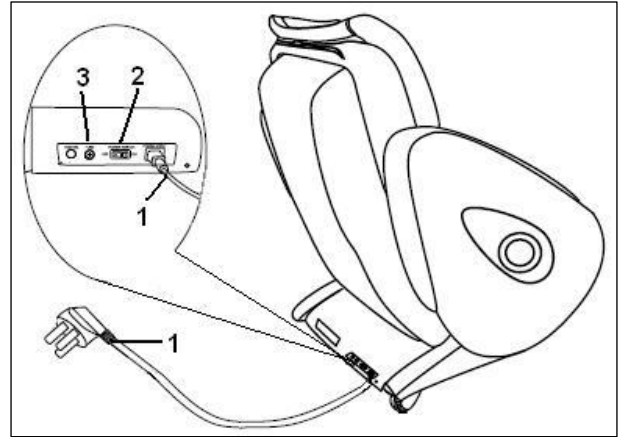Fig.17


## WARNING

Do not sit on the chair with the footrest raised.  
To avoid damaging the footrest, do not suddenly drop full body weight on to it.


### ❖ Reclining the backrest.

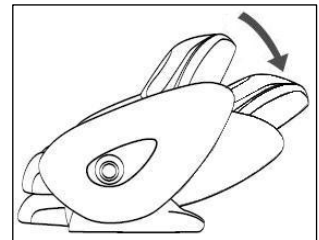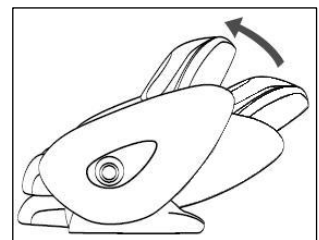- Press and hold the  key and the backrest will lower gradually.

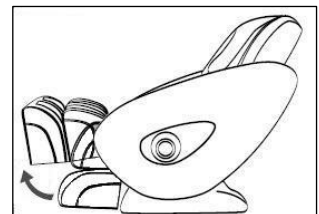
### ❖ Raising the backrest.

- Press and hold the  key and it will rise gradually.


### ❖ Raising the footrest.

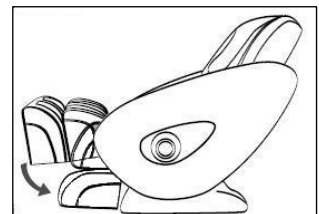- Press and hold the  key and the footrest will rise gradually.

### ❖ Lowering the footrest.

- Press and hold the  key and the footrest will lower gradually.

## ❖ Using the head massager.-


### 1. Assembly and installation

Place the massager on the backrest. Connect the quick-fit connectors C, D and E to those marked with the corresponding letters, Fig.18.

### 2. Methods of use:

- A. Find the best position for your massage.
- B. Lie back and rest your head.
- C. Press the massage button to start the massage.

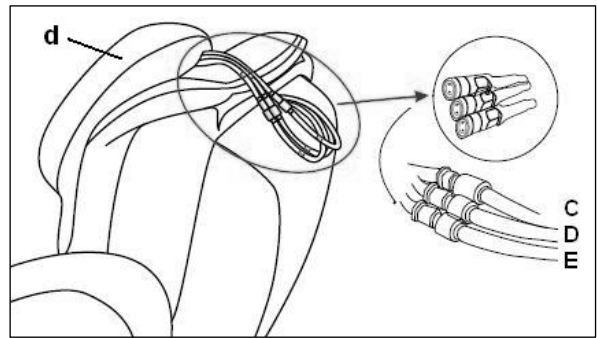Fig.18

## ❖ After use

- Once you have finished your relaxing massage, lower the legrest to the bottom position, Fig.19 and raise the backrest to its highest position, as shown in Fig.20.

Fig.19


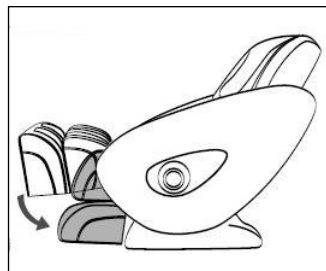
Fig.20
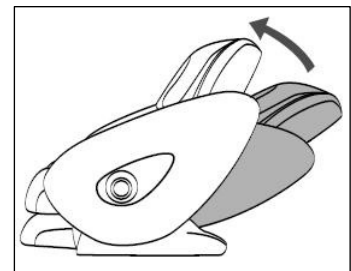
## ❖ After use.


- Press the ON/OFF  key and the chair will switch off. The massage rollers will return to their original position, set the switch located below the chair to the "O" position, Fig.21, and unplug the cable (1) from the wall socket.

- After use, switch the chair off and unplug it from the mains supply.

- Place the remote control unit in its support.

- To protect children from accidents and injuries, always make sure that the chair has been unplugged from the electrical mains socket.
- If the chair is not going to be used for prolonged periods, store it away.

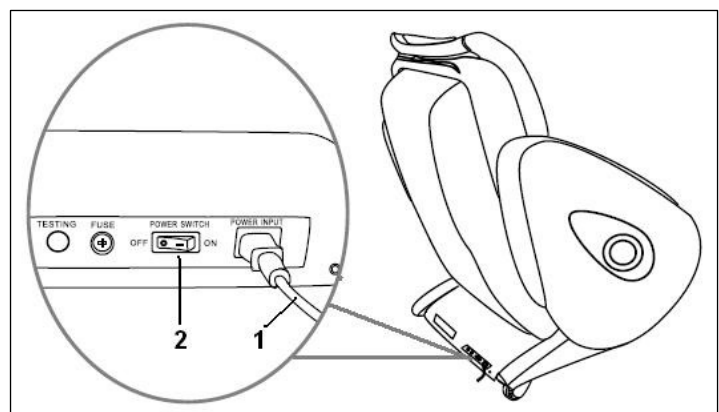Fig.21

# M1200

## ❖ FAULT FINDING

Problem	Check the following	Solution
Rasping sound during use	This is due to the construction of the chair or friction between the massage rollers and the leather upholstery,	Normal sound
It does not operate	Is it plugged into the electrical wall socket	Plug it in.
	It is switched off	Switch it on
	The fuse has blown	Replace the fuse with a similar type
The music does not play correctly	The music is not MP3 format	Change music format
The sound is echoed when you press the music key.	The disc has not been inserted correctly	Reinsert the disc
	The disc is still going through the start-up sequence	Wait 5-10 seconds before selecting any music
	Incompatible music format	The computer is checking the MP3 files
No sound from earphones	Make sure the earphones are plugged in.	Change the earphones.
	Volume is too low or switched off	Increase the volume
It is not possible to raise or lower the backrest or footrest	Is the backrest or footrest overloaded	Reduce the load or pressure on the backrest or footrest.
It is not possible to extend the footrest normally.	The footrest is too close to the floor, which is activating the floor proximity sensor.	Place the footrest in a higher position.

If the above solutions do not solve the problem, switch the chair off and unplug it from the wall socket, then get in touch with your Technical Assistance Service or local retailer, who will remedy the fault.

## • MAINTENANCE

TIP

Unplug the chair from the mains supply before performing any cleaning or maintenance tasks.

To prevent accidents do not attempt to unplug the chair with wet or damp hands or clothes.

Do not tug on the mains cable to unplug it.

### ❖ Maintenance of component parts.

#### Main body

Use a lightly dampened cloth with a mild cleaning product to clean the plastic elements and armrests  
Dry with a clean dry cloth.

#### WARNING

Do not use caustic cleaning products or detergents, etc. to avoid damage, fading, etc.

#### Remote control unit

Use a lightly dampened cloth with a mild cleaning product to clean the remote control unit

#### WARNING

To avoid accidents do not use moist products to clean the remote control unit.

#### Backrest cushion and pillow

Use a lightly dampened cloth with a mild cleaning product to clean these elements.  
Leave them to dry.

#### Storage

After cleaning, store the chair in a dry place. Cover it if it is not going to be used for prolonged periods.

#### WARNING

Do not expose the chair to direct sunlight or high temperatures as this will cause the surfaces to fade.

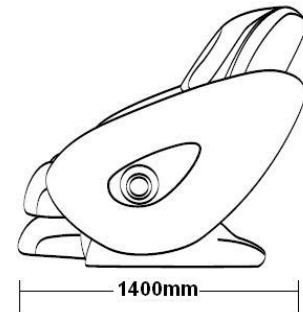
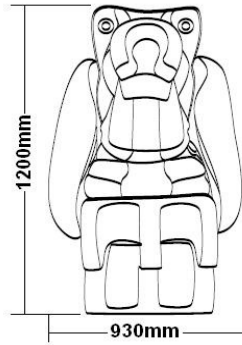
#### Inspection

 <b>Inspection</b>	<p><b>Following phenomena</b></p> <p>-</p> <ul style="list-style-type: none"> <li>-Burning smell</li> <li>-Hot electrical components</li> <li>-Other strange phenomena</li> </ul>	 <b>Stop use</b>	<p><b>In the event of accidents,</b> Unplug the chair and get in touch with the Technical Assistance Service</p> <p>–Do not dismantle any of the chair components or replace pieces that are not original manufacturers parts</p>
-----------------------	---	---------------------	---


## ❖ TECHNICAL DATA

- ❖ Product packaging: 1320x770x820mm.
- Side panel packaging: 960x610x380mm.
- Gross weight: 138kg.
- Net weight: 120kg.

- ❖ Dimensions in normal position
  - Length: 1,500mm.
  - Width: 930mm.
  - Height: 1,200mm.

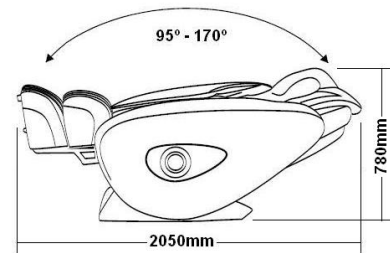- ❖ Extended dimensions
  - Length: 2,050mm.
  - Width: 930mm.
  - Height: 780 cm.

## ❖ Specifications:

- Model M-1200
- Voltage 220~240V AC
- Mains frequency 50~60Hz
- Power consumption 150W
- Timer 5/10/15/25/30 minutes
- Max. Time 30 minutes
- Maximum user weight: <120 kg
- Cable length:
  - Controller cable: 1.2 m
  - Mains cable: 2.1 m
- Conditions for use:
  - Room temperature: 20°- 60°
  - Humidity level: 30-85 HR
- Storage conditions:
  - Storage temperature: 20°-60°
  - Humidity level: 30-85 HR
- Safety function:
  - Safety: Class I
  - Equipped with overheat and overvoltage protection

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MAQUINASPORT, S.A.