**UNIDAD ELECTRÓNICA  
ELECTRONIC MONITOR  
UNITÉ ÉLECTRONIQUE  
ELEKTRONIKEINHEIT  
UNIDADE ELETRÓNICA  
IMPIANTO ELETTRONICO**

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
**902 170 258**

## ELECTRONIC MONITOR.-

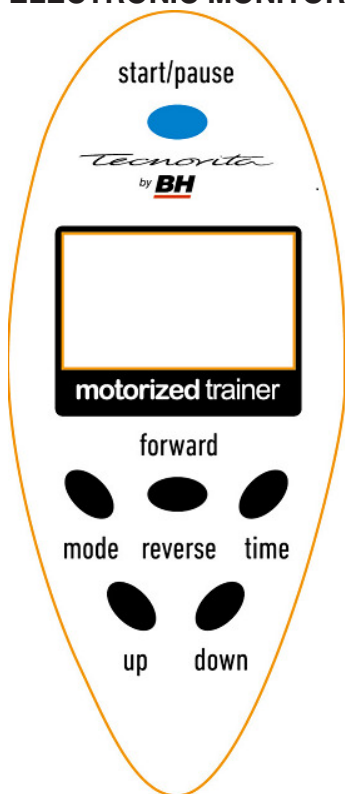**START/STOP**– On/Off. Press to start / finish the exercise.

**MODE**– Mode selection. Press the MODE key to select between P1, P2, P3, P4, P5 or manual mode. While you are exercising, press this key to display the time, distance, speed or all variables successively.

**FORWARD/REVERSE**– It changes the pedaling direction.

**TIME**– It controls the length of the exercise. Press it to increase the duration of the exercise.

**UP - SPEED (+)** – It increases the speed.

**DOWN - SPEED (-)** – It decreases the speed.

1. Plug the cord to the power supply and turn the switch located inside the motor box to the on position.

2. Press the **MODE** key to select the mode P1, P2; P3, P4, P5 or manual mode (-).

3. Press the **START / STOP** key to turn on the appliance.

4. After exercise press the **START / STOP** key.

5. Turn the switch located inside the motor box to the off position and and unplug the cord from the power supply.

### P1 MODE

Time 10 minutes, initial speed level 1, level change every 30 seconds (level 1 - 2 - 3 - 4 - 5). After 5 levels the pedaling direction and repetition speed rate change.

### P2 MODE

Time 10 minutes, initial speed level 2, level change every 15 seconds (level 2 - 3 - 4 - 5 - 6). After 5 levels the pedaling direction and repetition speed rate change.

### P3 MODE

Time 10 minutes, initial speed level 3, level change every 20 seconds (level 3 - 5 - 7 - 9 - 11). After 5 levels the pedaling direction and repetition speed rate change.

#### **P4 MODE**

Time 10 minutes, initial speed level 4, level change every 20 seconds (level 4 - 6 - 8 - 10 - 12). After 5 levels the pedaling direction and repetition speed rate change.

#### **P5 MODE**

Time 10 minutes, initial speed level 5, level change every 30 seconds (level 5 - 10 - 5 - 10). After 4 levels the pedaling direction and repetition speed rate change.

While MODES P1, P2, P3, P4, P5 are being used speed and pedaling dirección cannot be changed manually.

#### **MANUAL MODE (-)**

Time 5 minutes, adjustable pedaling direction, speed and time.

Do not hesitate to get touch with the Technical Assistance Service if you have any queries by phoning customer services (see last page in manual).

**BH RESERVES THE RIGHT TO MODIFY THE SPECIFICATIONS OF ITS PRODUCTS WITHOUT PRIOR NOTICE.**

## **Guarantee**

This product is guaranteed for 2 years in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25<sup>th</sup> May, 1999, regarding certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10<sup>th</sup> July; Guarantees on the sale of consumer goods.

To this effect, EXERCYCLE, S.L. VAT nº B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product.


Date of purchase.-	Purchased in.-
Purchaser.-	Dealer stamp.-

Please keep the receipt as proof of purchase.

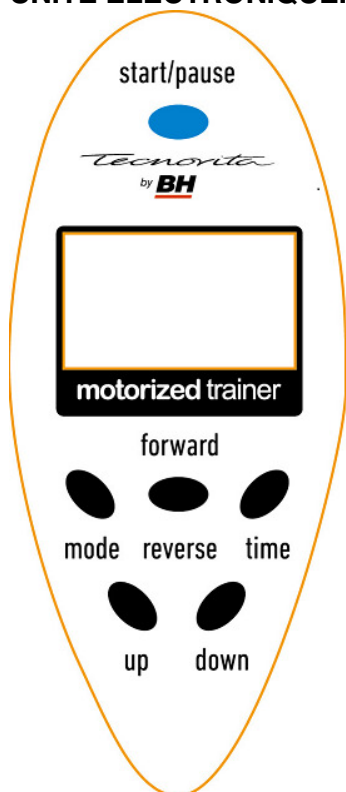Should you have any queries please do not hesitate to phone our customer service line:

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**BH Germany GmbH**